


Compact525

5 Litre Oxygen Concentrator

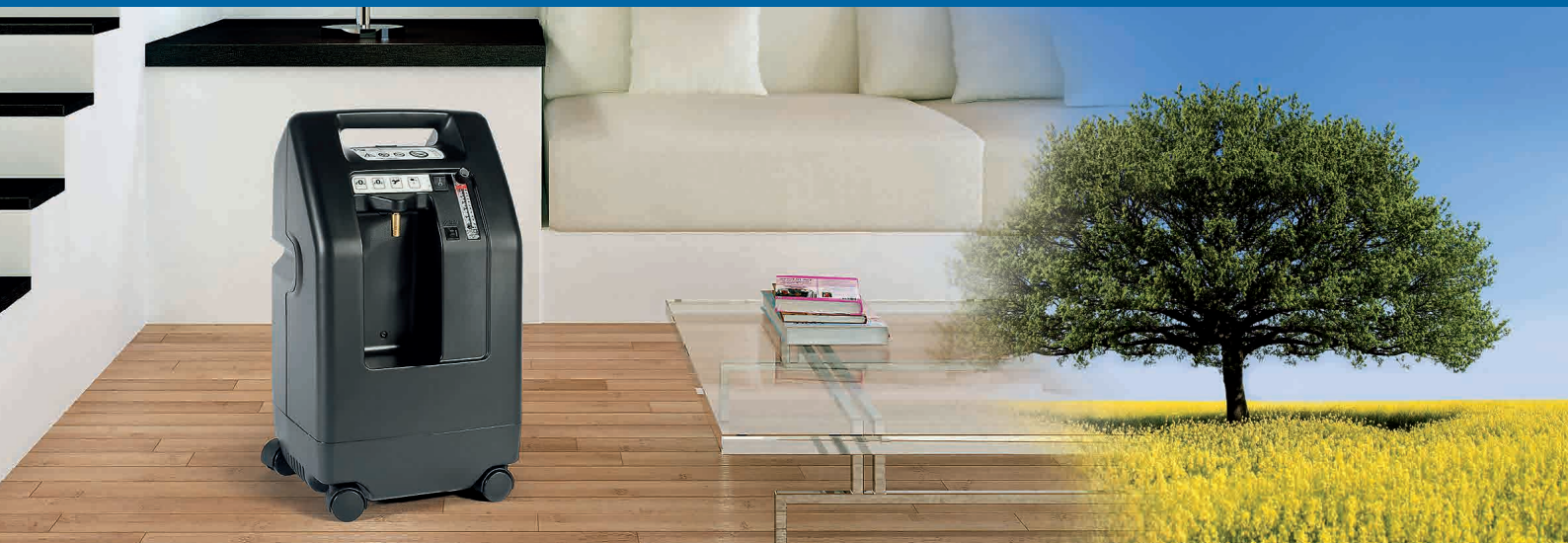

Compact525

Manufacturing products that are of the highest possible quality is key to DeVilbiss Healthcare. The latest Compact 525 oxygen concentrator is designed and built to meet the high demands of patients and cost challenges faced by service providers.

In commitment to our quality aims, we have invested in new manufacturing processes and innovative production techniques. The new black, Compact 525 is now manufactured at our dedicated oxygen therapy facility in Pennsylvania, USA, which has given us the control to deliver a reliable, high quality, more eco-friendly oxygen concentrator.

Delivering up to 5 litres per minute with high oxygen concentration across all

flow rates and adopting the latest developments in the reduction of power consumption, the latest Compact 525 increases patient comfort and lowers operating costs.

Added features

- Improved sound level reduction (525KS)
- Front label with easy to read pictograms
- New black colour

Service and maintenance friendly

- 2 part case design
- 3 years warranty as standard (an optional extended 2 year warranty package available)

Patient safety features

- OSD® Sensor for continuously monitoring oxygen levels
- Visual and audible alarms for low oxygen levels, power failure, pressure drop and service required
- Oxygen outlet incorporating a fire protection adapter

Environmentally friendly

- Intelligent power management system utilises Turn-Down Technology providing less power consumption below flow rates of 2.5 l/min
- Reduced packaging: reducing materials and waste, as well as transportation costs and emissions – without compromising product protection



SPECIFICATIONS

Dimension (H x W x D)

62.2 x 34.2 x 30.4 cm

Weight

16.3 kg

Flowrate

0.5 to 5 l / min

Oxygen Concentration

at 0.5 - 5 l / min

93% ± 3%

Electrical Requirements

230 V, 50/60 Hz

Power Consumption

approx. 290 Watt at 2 l/min.

approx. 312 Watt at 5 l/min.

Classification

II Type B

Models available

525KS:

Ultra-quiet 5 litre oxygen concentrator

525KS-LT:

5 litre oxygen concentrator



Compact 525



Front label with easy to read pictograms



Humidifier



Paediatric Flowmeter



Bar O2 Limited.
Unit 5A, Mount Industrial Estate,
Mount Road, Stone, Staffordshire,
ST15 8LL
United Kingdom
Phone: +44 1785 894000

Theoxygenstore.com

Distributor: Bar O2 Limited

Bar O2 Limited T/A the OxygenStore, Unit 5A Mount Industrial Estate, Mount Road, Stone, Staffordshire, ST15 8LL

Company Registered in UK No 04192727 VAT Registration No 770 4717 24
EORI GB770471724000 Dun and Bradstreet D-U-N-s No 232439047